

DENTAL WATERLINE SOLUTIONS BY HU-FRIEDY

Keeping waterlines clean is critical for every dental office. Evidence shows the narrow tubing, low flow rates and frequent stagnation periods found in water delivery systems create the ideal breeding ground for biofilm.¹

The CDC says dental practices should use water containing no more than 500 colony forming units (CFUs) per milliliter. However, with new waterlines, the microbial count can reach 200,000 CFUs/mL within just 5 days if left untreated*.

Hu-Friedy offers a variety of products that provide dental offices with an effective and economical way to both clean and maintain their waterlines.



Hu-Friedy Waterline Filter

The Hu-Friedy Waterline Filter continuously **works to deliver clean water to patients and staff for up to one year**. The filter features a convenient, ready-to-use cartridge that is easy to install and inhibits microbial contamination in the waterline. With a universal design, it can be used with either a municipal waterline or an independent water reservoir. The filter is non-corrosive and can easily be disposed of in the standard trash.

VistaTab™ Dental Waterline Cleaner Tablets

VistaTab Dental Waterline Cleaner Tablets are antimicrobial tablets used for the cleaning and control of microbial contaminants in dental unit waterlines. When placed in water, the solid tablet quickly dissolves to create a non-corrosive solution that is easy to use.

VistaClean™ Irrigant Solution Concentrate

VistaClean Irrigant Solution Concentrate is a daily solution derived from organic citrus botanicals and is safe for use with patients. It helps protect waterlines and equipment from deposits and scale and works to emulsify organic and inorganic contaminants that can build up in waterlines. VistaClean leaves no taste or odor and is non-toxic when used as directed. This **environmentally friendly solution is non-corrosive** and has no negative effect on bond strength.

* Reference: Centers for Disease Control and Prevention. Guidelines for Infection Control in Dental Health-Care Settings - 2003. MMWR 2003; 52 [No. RR-17]: [p.28]

How the best perform



THE CDC RECOMMENDS...

that dental unit waterlines be tested on a regular basis to maintain its established standard of less than 500 colony-forming units per milliliter.*

If dental unit waterlines are not maintained on a regular basis, bacteria will develop on the inner surfaces of the tubing, blocks and other system components.

Simply maintaining the waterlines is not enough; it is a two part process of cleaning and maintenance. Over time, bacteria can overwhelm the waterline environment and begin to create microbial deposits. Both maintaining with a product like Hu-Friedy Waterline Filter or VistaClean, along with cleaning the waterlines regularly with an antimicrobial, like VistaTab, will clean and control microbial contaminants to provide improved water quality while helping to protect patients, staff and equipment.

* Reference: Centers for Disease Control and Prevention. Guidelines for Infection Control in Dental Health-Care Settings - 2003. MMWR 2003; 52 [No. RR-17]: [p.31].

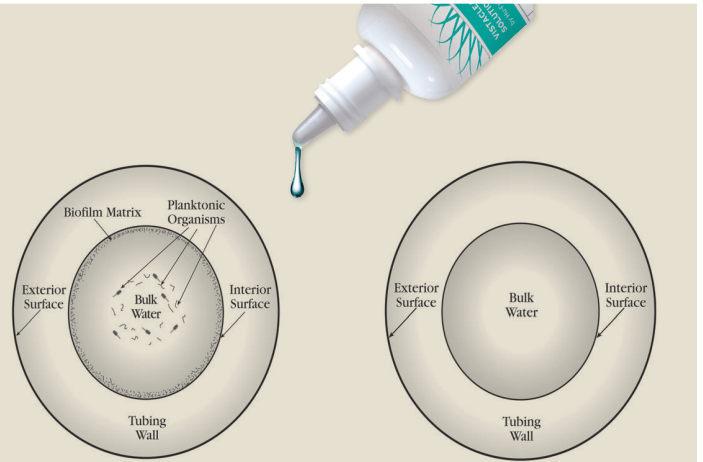


Figure 1. Contaminated dental waterline. Biofilm on interior surface with abundant free-floating organisms.

Figure 2. Clean dental waterline tubing. No biofilms or free-floating organisms present.



DENTAL WATERLINE PRODUCTS



Hu-Friedy Waterline Filter[†] IMS-1440

The Hu-Friedy Waterline Filter is an annual dental unit waterline treatment filter that inhibits microbial contamination and odor causing bacteria.

VistaTab is EPA registered in all 50 states. VistaClean is an FDA cleared irrigant.
†Available only in the U.S.



VistaTab Dental Waterline Cleaner Tablets[†] IMS-1451

This box contains: 15 Tablets (1.7 grams each)
Net Weight: 22.5 grams (0.90 oz.)



VistaClean Irrigant Solution Concentrate IMS-1452

One Fluid Ounce Bottle (29.57 ml)

Visit us online at www.Hu-Friedy.com

1. Center for Disease Control and Prevention, 2018. Dental Unit Water Quality. <https://www.cdc.gov/oralhealth/infectioncontrol/summary-infection-prevention-practices/dental-unit-water-quality.html>

Call 1-800-Hu-Friedy or contact your authorized Hu-Friedy representative for more information.

©2020 Hu-Friedy Mfg. Co., LLC. All rights reserved. All company and product names are trademarks of Hu-Friedy Mfg. Co., LLC its affiliates or related companies, unless otherwise noted. IMS-337/0520

